

# Xcela-250 Series™

## CO<sub>2</sub> PURITY ANALYZER

PRODUCT BRIEF  
TURN-KEY TABLE RACK CO<sub>2</sub> ANALYZER



### KEY BENEFITS

- Comprised of rugged, reliable analyzers with totally Silconert® passivated hardware & tubing that meet all ISBT guidelines for sampling & testing of beverage-grade CO<sub>2</sub>.
- This Xcela™ system offers *on-line* measurement of Total Sulfur Content (TSC), Aromatic Hydrocarbon (Benzene), Acetaldehyde (AA), Methanol (MeOH), Total Hydrocarbon Content (THC) and Water Vapor (H<sub>2</sub>O).
- ALI's Methods Manual includes all validated ISBT / CGA / EIGA specified methods. All systems are shipped fully plumbed, tested & professionally crated for safe transport plus rapid on-site start-up.
- An ASCM-10 gas controller allows for either operator or automated data station control of up to 10 sample lines plus analyzer auto-verify / auto-cal functions.
- An intuitive, industrial-grade data station (DSA Blackbox™) conveniently displays up to 32 channels of test data, performs data-logging, allows manual test data input (ex. DT, NVR, % CO<sub>2</sub>, Sensory), issues quality alarm alerts, generates CoA reports (for several purity grades), features operator-programmed sample-line sequencing & remote access for data quality auditing or preventive maintenance services.
- ALI, along with its worldwide service network offer installation, training, expert technical support, ISO-certified lab back-up, preventative maintenance, repairs, spares, all calibration & support gases plus all consumables (ex. all Mfg DT models).



**BASILINE BAM 8900**  
Aromatic Hydrocarbon (Benzene)  
Acetaldehyde (AA), Methanol (MeOH)



**BASILINE TSC 8900**  
Total Sulfur Content Analyzer



**GE DEW IQ**  
Moisture Analyzer



**ASCM-10-2**  
Automated Sample Line and  
Calibration Control Module



**BASILINE THA 9000**  
Total Hydrocarbon Analyzer (THA)



**DSA BLACKBOX™**  
Data Station

### RECOMMENDED ACCESSORIES




Hydrogen Generator  
Proton On-Site  
Hogen-600



Zero Air Generator  
Parker Balston  
with GAST Compressor  
& Auto-Drain.



PowerVar  
Power Conditioner &  
Surge Protector

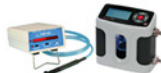

Digital Flow Meter  
(Pressure &  
Temperature Compensated)  
& Mini Leak Detector



Calibration Gases  
(available in various sizes)



Detector  
Tubes



Sampling Cylinder  
with Vaporizing Pressure  
Regulator & stand.  
(for Liquid CO<sub>2</sub> Sampling)



ZCheck 5XLT  
Passivated Check Valve Assy  
& Zahm & Nagel % Purity



**Airborne Labs International**

ISO / IEC 17025:2005

**AIRBORNE LABS INTERNATIONAL, INC**

22C World's Fair Drive, Somerset, New Jersey 08873  
T: (+1) 732-302-1950 F: (+1) 732-302-3035  
E: AirborneLabs@aol.com www.AirborneLabs.com